**emmi**<sup>®</sup> ULTRASONIC

**emmi**<sup>®</sup> DENT **emmi**<sup>®</sup> PET

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**READ  
CAREFULLY!**

**QUICKGUIDE +  
PROBLEM FIXING**

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The device is working correctly.

**Red LED:**

The device needs to be charged or the attachment is defective.

Replace attachments **every 3 months**.


### Video tutorial:

Correct placement of the ultrasound attachments for optimal ultrasound transmission!

Scan me
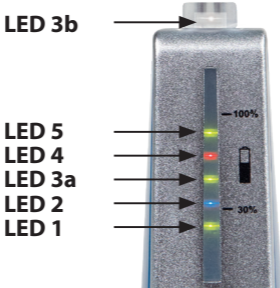
### Functions of the LEDs:

#### a) Charging control function LED2

**Blue LED:** emmi®-dent is on the charging station

**Light blue light:** emmi®-dent in charging process

**Weak blue light:** emmi®-dent is fully charged

#### b) Battery status

##### (LED1, LED3a, LED5 & LED3b)

The three **green** LEDs (and the white LED on the brush head) are used in the manner of a sequential light during the cleaning process.

Is the battery fully charged, the LEDs light up as follows:

- LED1, LED3a, LED5 + LED3b

If the battery is approx. 1/3 discharged, the following LEDs light up in sequence:

- LED1, LED3a + LED3b

If the battery is approx. 2/3 discharged, only the following LED lights up:

- LED1

#### c) „Urgent recharge“ indicator

If the battery is heavily discharged, the LED4 flashes **red** a few times after the cleaning process. Please charge the emmi®-dent.

#### d) Switching between vibration and acoustic signal (except emmi®-dent Platinum Care)

If the On button of the emmi®-dent is pressed for approx. **10 seconds**, the LED2 flashes **red** and the vibration as well as the acoustic signal are switched off.

Press and hold down the On button for approx. **10 seconds**, and the LED2 flashes **green**. The vibration as well as the acoustic signal are switched on now.

Before first use charge for **24 hours**.

Insert the attachment **vertically**.

**Green LEDs:**

The device is working correctly.

**Red LED:**

The device needs to be charged or the attachment is defective.

Replace attachments **every 3 months**.


**Video tutorial:**

Correct placement of the ultrasound attachments for optimal ultrasound transmission!

Scan me
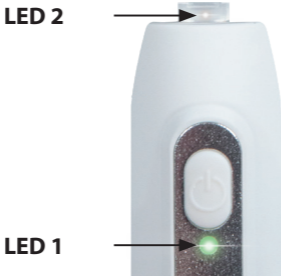
### Functions of the LEDs:

#### a) Charging control function LED2

**Blue LED:** emmi®-dent is on the charging station

**Light blue light:** emmi®-dent in charging process

**Weak blue light:** emmi®-dent is fully charged

#### b) Battery status (LED1 & LED2)

The two LEDs are used as battery status indicators during the cleaning process.

If the battery is between full and approx. 1/3 discharged:

- LED1 lights up permanently **green** and LED2 once approx. every 3 seconds

If the battery is approx. 2/3 discharged, only the following LED lights up:

- LED2 lights permanently **green**

#### c) „Urgent recharge“ indicator

If the battery is heavily discharged, the LED1 flashes **red** a few times after the cleaning process. Please charge the emmi®-dent.

#### d) Switching between vibration and acoustic signal (except emmi®-pet)

If the On button of the emmi®-dent is pressed for approx. **10 seconds**, the LED2 flashes **red** and the vibration as well as the acoustic signal are switched off.

Press and hold down the On button for approx. **10 seconds**, and the LED2 flashes **green**. The vibration as well as the acoustic signal are switched on now.

#### Please consider:

emmi®-pet is completely soundless and vibrationless.

Here you will find independently correctable errors that could occur with your emmi®-device:

### **1) Flashing of the red LED means:**

#### **a) Battery discharged:**

Place the handpiece on the charging station so the emmi®-ultrasound device can start the charging process.

#### **b) Indication of a defective attachment:**

Place a new attachment on the handpiece and switch the device on again. The **green LED** should now light up continuously.

### **2) Green LED lights up and turns off immediately:**

**a) No attachment attached:** Please attach an attachment to the handpiece and switch it on again. The **green LED** should now light up continuously.

**b) Attachment defective:** Please attach a new attachment to the emmi®-ultrasonic device and try to switch it on afterwards. The **green LED** should now light up continuously.

### **3) The green LED lights up continuously, but without vibration:**

#### **Vibration was manually switched off:**

Please place an attachment on your emmi®-device. Press and hold the **ON/OFF key** for approx.

**10 seconds** to switch on. The LED lights up **green three times**. You can now remove your finger from the ON/OFF key. After switching the device on, the vibration should work again.

To **switch off** the vibration, please proceed in the same way. The LED will light up **red three times**.

(Please note that our models emmi®-dent Platinum Care and emmi®-pet do not have a vibration motor and therefore this function cannot be used with these models).

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**E-Mail:** [info@emmi-ultrasonic.de](mailto:info@emmi-ultrasonic.de)

### **You need personal assistance?**

Please feel free to contact our customer service:

**Phone:** +49 (0) 6105 / 40 67 94

**WhatsApp:** +49 (0) 160 / 50 91 446

**E-Mail:** [info@emmi-ultrasonic.de](mailto:info@emmi-ultrasonic.de)

Emmi Ultrasonic GmbH

Gerauer Str. 34 | 64546 Moerfelden-Walldorf

Tel.: +49 (0) 6105 / 40 67 94 | [info@emmi-ultrasonic.de](mailto:info@emmi-ultrasonic.de)