emmi[®]DENT

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emmi[®]-dent Platinum

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Emmi Ultrasonic GmbH • Gerauer Strasse 34
64546 Mörfelden-Walldorf • Germany
Tel: +49 (0)6105 - 40 67 94
Email: info@emmi-dent.de

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emmi®-dent - The original device for innovative dental care with 100% Original Ultrasound

We thank you for your purchase of an emmi®-dent Platinum ultrasound device for your dental cleaning and oral hygiene („ultrasound toothbrush“). Please read carefully this user manual for optimal results.

Content emmi®-dent Platinum Set

1 emmi®- dent Platinum hand piece,
1 emmi®- dent ultrasound attachment for adults,
1 emmi®- dent charging station,
1 emmi®- dent tube ultrasonic toothpaste (75ml)

Loading and service life per load

Important – Charging usage term

For each charging, connect the plug of the charging station into the socket, insert the hand-piece into the dock. When the charging station is connected properly, the green LED light on the charging station is ON. Charge the hand-piece for 24 hours for the first time. **Further information can be found in the LED display section (p.13).** If you use the device 6 minutes daily, you can use it with the battery for approx. 15 days. Charge it then as needed.

Before the first application

Insert the ultrasound brush head vertically in the hand-piece. Also remove the ultrasound brush head in vertical direction.

Switching on and off

Switching on:

Switch on the device by pressing the key (for approx. 1 second). To switch off the vibration motor and the signal tone, please press and hold the On/Off key for approx. 9 seconds. The LED lights up red four times. To switch on the vibration motor and the signal gate again, please keep the On/Off key pressed again for approx. 9 seconds. The LED lights up green four times.

The green LED of the power indicator lights up. The battery charge level is indicated by two more green flashing LEDs. Furthermore, a vibration can be felt, which represents an additional function message for the user.

Switch off/Automatic switch off:

Switch off the unit by pressing the push-button (for approx. 1 second). The unit switches off automatically after approx. 6 minutes (Battery Save Function). The device is also provided with an electronic sensor. If a wrong or defective ultrasound brush head is used, or no brush head has been inserted, the unit will also switch off automatically, the green LED goes off and the red LED blinks. Don't insert the brush head during transport, the device can not be accidentally switched on (Travel Save Function).

LED Display

Functions of the LEDs:

a) Charging control function LED2

Blue LED: emmi[®]-dent is on the charging station

Light blue light: emmi[®]-dent in charging process

Weak blue light: emmi[®]-dent is fully charged

b) Battery status

(LED1, LED3a, LED5 & LED3b)

The **three green LEDs** (and the white LED on the brush head) are used in the manner of a sequential light during the cleaning process.

LED 3b

(from mid 2022)


LED 5

LED 4

LED 3a

LED 2

LED 1

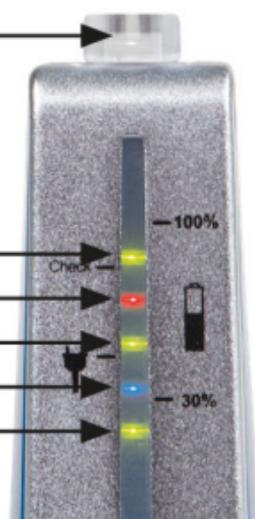Is the battery **fully charged**, the LEDs light up as follows:

• LED1, LED3a, LED5 + LED3b

If the battery is approx. **1/3 discharged**, the following LEDs light up in sequence: • LED1, LED3a + LED3b

If the battery is approx. **2/3 discharged**, only the following LED lights up: • LED1

c) „Urgent recharge“ indicator

If the battery is heavily discharged, the LED4 flashes **red** a few times after the cleaning process. Please charge the emmi[®]-dent.

d) Switching between vibration and acoustic signal (except emmi[®]-dent Platinum Care)

If the On button of the emmi[®]-dent is pressed for approx. **10 seconds**, the LED2 flashes **red** and the vibration as well as the acoustic signal are switched off.

Press and hold down the On button for approx. **10 seconds**, and the LED2 flashes **green**. The vibration as well as the acoustic signal are switched on now.

Using the emmi[®]-dent Platinum with 100% Original emmi[®]-dent Ultrasound.

After switching on the unit you will feel a soft vibration and hear a gentle sound. The actual cleansing process is not affected by this. You cannot see, hear or feel ultrasound and thus a small motor is integrated within the hand piece, so that the user can feel and hear that the device is switched on and ready for operation. Also, a low peep sound is emitted every 6 seconds and the motor stops briefly to give you a haptic feedback. (Beep can be switched off manually).

Daily use

For your dental hygiene, we recommend you to clean your teeth with the 100% Original ultrasound of your emmi®-dent Platinum twice a day for approx. 3 min.

1. Moisten the bristles of the ultrasound brush shortly under running water.
2. Apply a double pea-sized quantity of emmi®-dent ultrasound toothpaste on the bristles of your ultrasound brush head. Only the emmi®-dent ultrasound toothpaste guarantees the Original 100% emmi®-dent ultrasound function. For more details on the emmi®-dent ultrasound toothpaste, please read under the separate section, "emmi®-dent ultrasound toothpaste" in this user manual.
3. Switch on the device and bray the ultrasound toothpaste uniformly on the inside and outside surfaces as well as on the chewing sides of your teeth.
4. The bristles of the ultrasound brush head cover 2 – 3 teeth, according to the size of your teeth. Maintain the bristles gently and without pressure motionless on the inside, outside and on top (chewing sides) of the teeth for approx. 6 seconds each time. The actual duration depends on the amount of impurities on the teeth (plaque, nicotine, wine stains, etc.). In case of acute problems, such as gingivitis or parodontitis, leave the bristles for approx. 15 to 20 seconds on the affected sites. This will result in a quicker and better healing. A longer application does not harm your teeth, because the emmi®-dent Platinum implements a gentle ultrasound. When you have finished your dental hygiene, move the bristles various times over your teeth, in order to remove any further adherences or impurities. Then rinse your mouth with water or an oral rinse.
5. After cleansing, wash the bristles and hand-piece under running water. Switch off the device and dry the hand piece. Please note that the bristles are free and clean from leftover food particles or other impurities. Now you can place the device with the integrated ultrasound brush head back on the charging dock.

1.

2.

3.

4.

5.

Inflammations, Parodontitis and Mouth Ulcers

In acute cases, it is recommended to hold the bristles gently for 30 to 60 seconds on the affected sites. This will result in a quicker and better healing. The Original emmi®-dent ultrasound is harmless and can also be used for a longer period of time.

Bad Breathe

Bad breath (halitosis) is caused by living bacteria in the mouth. The emmi®-dent Platinum Original ultrasound breaks down the bacteria on and in the gums and bad breath disappears for up to 12 hours.

emmi®-dent ultrasonic toothpaste

emmi®-dent ultrasound toothpaste is available in five versions: emmi®-dent Whitening, emmi®-dent Fresh, emmi®-dent Kids, emmi®-dent Nature and emmi®-dent Nature Calcium. The both last mentioned ultrasound toothpastes is fluoride-free and vegan. All emmi®-dent ultrasound toothpastes are gluten-free and do not contain any nanomaterials.

It is often thought that ultrasound would require water as a transmission medium, but this is a false concept. Actually, the correct transmission medium depends on ultrasound frequency and ultrasound performance. It is a well-known concept in medical implementations. In this domain, creams and gels are used as adequate transmission medium. For your dental hygiene with Original emmi®-dent ultrasound, you use a special toothpaste (emmi®-dent ultrasound toothpaste). This toothpaste stands out for its ultrasound feature, because emmi®-dent ultrasound toothpaste is the only medium in which millions of micro-bubbles are located. Conventional toothpastes in the market only generate a little amount of micro-bubbles. Always consider an emmi®-dent ultrasound toothpaste for your dental care with the emmi®-dent ultrasound devices.

Replacement ultrasonic attachments

The following ultrasonic attachments are available:

emmi[®]-dent M2/M2 Wave/M4 (adults),

emmi[®]-dent S2/S4 (children or women),

emmi[®]-dent PM2 (clasp carrier).

Further information is available at

www.emmi-dent.com/equipment/

Please take into account that according to dentists, the ultrasound brush head should not be used by more than one person. One and the same ultrasound brush head must always be used by one and the same person. When changing the ultrasound brush head, ensure that the bristles of the brush head always point forwards in the direction of the ON/OFF button. Spare brush heads can be purchased at the locations where you have bought your emmi[®]-dent Platinum.

How it works

Our 100% Original ultrasound devices are the most innovative and effective method for a perfect dental care and oral hygiene (nano-cleansing). The emmi[®]-dent ultrasound devices have been developed by German engineers according to the latest scientific findings. Other than customary toothbrushes with mechanic friction and abrasive toothpaste substances (macro-tooth cleansing), they set the trend for a new way of motionless, gentle cleaning without brushing or friction. For decades, ultrasound technology has been implemented with great success for medical purposes. It has been applied for brain, oftalmologic, pregnancy and stomach research.

Ultrasound for dental and mouth hygiene is easy to use, whether we are young or old, because it is absolutely harmless and innocuous. Also carriers of pacemakers can use it without any risk. Research has proven that ultrasound of up to 50 watt is harmless for the human body. The emmi[®]-dent Platinum works with only 0.2 watt. The performance and safety of ultrasound waves have been tested by dentists and implantologists. Many years of experience have never reported any negative effects, also for crowns, inlays, implants or veneers. The ultrasound waves are created inside an ultrasonic microchip (Piezo Chip) which is patented worldwide. This chip is embedded inside the brush head and transmits the waves via the bristles to the ultrasound toothpaste medium. Millions of micro-bubbles are generated. Due to their microscopic size they are able to reach into the smallest spaces, gaps and fissures between the teeth and into the gum pockets. This also happens for wearers of implants and dental braces. The Original ultrasound of the emmi[®]-dent Platinum for gentle dental cleaning and oral hygiene works without any mechanical pressure, friction or strain. Due to the powerful implosion of the micro-bubbles generated by the air oscillations of the ultrasound, all impurities such as plaque or leftover food particles as well as bacteria and germs are removed efficiently. Ultrasound waves have also a positive effect on gums, gingiva and the tongue.

emmi[®]-dent Platinum cleans where your toothbrush or floss can't reach. The ultrasound removes bacteria and germs deep in the gums.

The positive effect of Ultrasound releases also body's own healing energies, it facilitates lymphatic drainage and activates the blood flow because emmi®-dent Platinum Original ultrasound penetrates deeply into the gum pockets. The prophylactic and healing effect of dental cleansing and oral hygiene with the Original emmi®-dent ultrasound is generated by the effective removal of most germs and bacteria. Not only for teeth and gums, but also inside the gum tissue.

Battery-charging technology

After use, the hand piece can be placed back in the charging dock for recharging the lithium battery or for a full battery charge for 15 days of toothbrush pleasure. If you use your emmi®-dent Platinum for more than 6 minutes per day, e.g. in case of several users, take a shorter interval for recharging the battery. If the battery is nearly empty, the power indicator does not light up when switching on the emmi®-dent Platinum (bottom green LED). In this case the red LED shortly flashes for 5 seconds and the device is switched off completely. After recharging the device for 24 hours, the normal functions are restored.

Battery exchange and disposal

The unit contains a special, environment-friendly, long-lasting battery. The battery is contained in the handpiece of the device. If after long-term use of the unit, the red LED flashes in spite of extensive recharging, the battery must be exchanged. Only use an emmi®-dent Platinum Lithium-Save-battery for the emmi®-dent Platinum (original exchange part with internal safety switch). If you cannot exchange the battery yourself, we will be happy to do this for you in a professional manner with low costs. If you exchange the battery yourself, please do the following: Loosen the screw at the bottom of the handpiece and remove the battery holder jointly with the bottom lid. When inserting a new battery, check the correct polarity. Insert the positive terminal in the bottom lid holder. During assembly, please ensure that the parts are correctly seated and sealed. Only then, safety class IP54 is guaranteed. You can obtain replacement batteries exclusively from Emmi Ultrasonic GmbH.

Note: The warranty, if still applicable, is only valid when the product is handled properly and the original spare parts are used.

The special battery in the device may not be disposed as normal domestic waste. For disposal, community collection points as well as other containers are available (e.g. at all sale points for similar batteries). Please see the hazard instructions mentioned below.

Hazard instruction

The emmi®-dent Platinum contains a battery with build-in safety technology. Even then, the battery may not come to kids' hands and the instructions for disposal must be followed. If in spite of all safety measures batteries are leaking, avoid contact with your skin and dispose the battery as described above. Safety measures are also required for the connection to the power supply (socket). The charging dock has a fully insulated cast housing, but never immerse the charging dock in water. The connecting cable must be in flawless condition. In case of damage to the housing, plug or connections, the whole charging dock must be discarded and replaced. Only use the emmi®-dent Platinum for its intended use as described. Children may only use the device according to the instructions and with the supervision of an adult.

Warranty

The emmi®-dent Platinum has a warranty of 24 months from the date of purchase, and a warranty of 6 months for the battery from the date of purchase. Any defects in the material and/or manufacturing errors that occur during this warranty term will be resolved without any costs. Your unit will be repaired or replaced. To the extent permitted by law, all other claims are excluded. The warranty will not apply in case of interference or dismantling by third parties without the prior approval of

Emmi Ultrasonic GmbH. The warranty will also be void for products that have been damaged through improper handling or incorrect use. If under a warranty repair is granted, the guarantee term will not be extended. For claims that are not justified, for example, in case of installation or operating faults, we reserve the right to claim the costs incurred for the same.

For claims under the warranty it is essential to provide the sales receipt stating the date of the purchase. For the fast processing of your claim, please send the complete unit including brush head and charging dock with sufficient shipping payment to Emmi Ultrasonic GmbH. Unfortunately, unfree shipments cannot be accepted.

Herewith we point out the legal rights of the consumer as well as that these are not limited by the warranty. We thank you for your understanding.

Intended use

The products may only be used as described in the user manuals and technical descriptions, which is for dental hygiene. Any use of non-original parts or components is subject to the recommendations or approval by Emmi Ultrasonic GmbH. Any other use or which goes beyond that use is considered as improper use. Emmi Ultrasonic GmbH accepts no liability for any damages that may result from the non-intended use of the product, and we assume no warranty for the error-free and correct functioning of the components. Intended use also includes the observation of all instructions in the user manual. The symbol on the type plate indicates the need to consult the operating instructions provided with the product. Trouble-free and safe operation of the product can only be ensured if the measures relating to proper transport, storage, assembly, installation and careful operation and maintenance are strictly observed. The tests according to BGV A3 must be performed according to the prescribed intervals. During the installation of the unit, ensure free access to the power socket on the charging dock.

Failure to use the product according to its intended use

Improper use of the product may result in hazard for life and limb as well as material damage. For this reason any use of the product other than according to its intended use is strictly prohibited. Emmi Ultrasonic GmbH accepts no liability for any damages that may result from improper use and assumes no warranty for the error-free and correct functioning of the components. For the proper operation of the unit, the following instructions must be followed:

- Only use the emmi®-dent Platinum for your dental cleaning and oral hygiene with the adequate voltage as described on the type plate (at the bottom of the charging dock)
- Do not use flammable or aggressive chemical cleansing products or acid solutions, etc.
- If damages due to transport are detected, do not connect the charging dock to the power supply.
- Only for indoor use at ambient temperature.

© Copyright EMAG AG. All rights reserved.
Manufactured in Germany/EU. The contents of these operating instructions are subject to copyright. The operating instructions may not be copied in whole or in part. Modifications, abridgements, extensions and supplements as well as any publication or translation by third parties require the prior written consent of EMAG AG.

emmi® ULTRASONIC

Emmi Ultrasonic GmbH • Gerauer Strasse 34
64546 Mörfelden-Walldorf • Germany
Tel. +49 (0)6105 40 67 94 • Fax +49 (0)6105 40 67 81
www.emmi-dent.de • info@emmi-dent.de

RoHS

(§6Abs.3, ElektroG)

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© Copyright EMAG AG. Alle Rechte vorbehalten.